


# Magnetic Resonance Imaging



Mark Cohen  
 UCLA Center for Cognitive Neuroscience  
 Departments of Psychiatry, Neurology, Radiology, Psychology,  
 Biomedical Physics, Biomedical Engineering

## The Plan


- The Magnetic Resonance Phenomenon & Contrast
- Spatial Encoding
- The "Pulse Sequence" Rules Everything  
*Seventh Inning Stretch*
- Functional MRI
- Image Quality and Artifacts



## Metaphor


['metəfɔːr, -fər]  
 noun

1. a figure of speech in which a word or phrase is applied to an object or action to which it is not literally applicable: *"I had fallen through a trapdoor of depression," said Mark, who was fond of theatrical metaphors | her poetry depends on suggestion and metaphor.*
  - a thing regarded as representative or symbolic of something else, esp. something abstract: *the amounts of money being lost by the company were enough to make it a metaphor for an industry that was teetering.*
2. Little white lie





## Spin

- Similar to Angular Momentum
- Intrinsic Property of Matter (rotation *not* required)
- Spin takes half-integer (0, 1/2, 1, 1-1/2) quantum values
  - Particles with integer spin are *Bosons*
  - Particles with half integer spins are *Fermions*
- Atoms with odd Mass Numbers (<sup>1</sup>H, <sup>13</sup>C, <sup>19</sup>F, <sup>23</sup>Na) have half-integer spin.
- Electrons, Neutrons, Protons have spin = 1/2





## Spin States

- A Spin 1/2 particle has two states ("up/down", "1 and 2",  $\alpha/\beta$ )
- In a magnetic field,  $B_{\text{ext}}$ , the two states have different energies

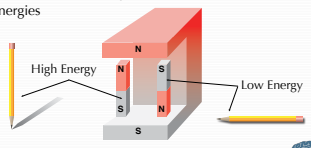



## Transition to Equilibrium

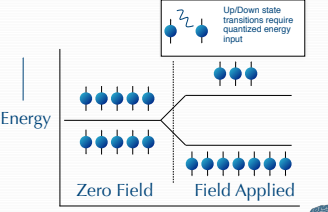When the nuclear magnet is placed into a magnetic field it can adopt only one of two states - aligned with or aligned against - the external field. Then the individual nuclear spins are aligned AGAINST the external field they are in a lower energy state. Moving a system to its lowest energy state is called "relaxation" and amounts to a loss of stored heat energy.

Absent an external field, there will be an equal number of spins in the two allowed states. When an external field is applied, some of the spins change state to oppose it. The local effects on individual spins is much larger than the effect of an external field; therefore only a small number of spins change state (1 in a million or so). The state changes actually require the capture of energy from the environment, which occurs rarely. Therefore, the relaxation to magnetic equilibrium can take some time.

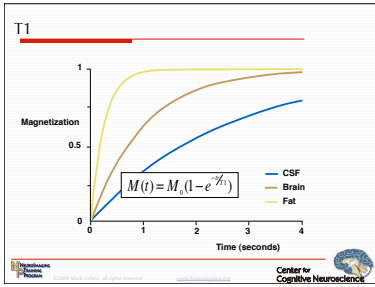

The approach to magnetic equilibrium is exponential. When the system is far from equilibrium, most interactions with the environment will move individual spins to the lower energy state. As the system reaches equilibrium, however, the state transitions become nearly equally likely to be either up or down. The exponential process is governed by a single time constant, T1, which is characteristic of the tissue type. The T1 of water (CSF) is up to several seconds, while the T1 of lipid (fat) is a few tenths of seconds. Most body tissues lie between. The difference in T1 forms a major contrast mechanism in MRI.

### Nuclear Spin

- Angular momentum is a vector quantity having magnitude & direction.
- Angular Momentum is Conserved, causing precession to occur

Precession Axis  
Spin Axis

Center for Cognitive Neurosciences

Nuclear spin is a quantized property. The quantum angular momentum behaves like the classical spin of a moving top. Unpaired nuclear spin results in the creation of a magnetic moment, oriented along the axis of spin.

### Proton Precession

Due to their angular momentum, Protons precess in the magnetic field.

Applied Magnetic Field: B

Precession:  $\omega$

Spin

$\omega = \gamma \times B$   
 $\gamma_H = 267.52 \text{ Rad/sec/Tesla}$   
 $= 42.577 \text{ MHz/Tesla}$

Center for Cognitive Neurosciences

The angular moment of the spins is conserved. This means that when an external field is applied, the individual spins will precess about it. The rate of precession is known as the Larmor frequency and is proportional to the strength of the external field. The "Larmor constant" or "gyromagnetic ratio" is about 42.577 million cycles/second per Tesla.

### MRI 101a

Protons, the Nucleus of Hydrogen, Have a Magnetic Moment

MAGNETIC MOMENT

SPIN

Center for Cognitive Neurosciences

To understand how the ultra-low field imager can work, it is useful to review a very few basics of MRI. MRI works by detecting the small magnetic moments of atomic nuclei, principally protons, which are in high abundance in the human body.

### MRI 101b

Protons Align (polarize) and Precess in Magnetic Fields

Center for Cognitive Neurosciences

When placed in a magnetic field, the protons align with the applied field and "precess" about it. The strength of the alignment, which is called the polarization, determines the strength of the MRI signal that will be received. This polarization shows up as a small excess of protons aligned along the field, rather than aligned against it.

The precession of the protons creates a time-varying magnetic field, that is detected by the imager.



When an ensemble of protons is placed into a magnetic field each of the individual spins adopts the spin up or spin down condition. The precession of the spins also implies that each has a rotational phase. In general, the phases of the individual spins can be expected to be random with respect to one another.

Each spin has a Longitudinal vector component about which it precesses and a Transverse component that rotates about the applied field.

### Proton Responses to Magnetic Field


- Spin Alignment Along Net Applied Field  
spins align parallel or anti-parallel to the applied field
- Precession About the Magnetic Field  
at a precession frequency of:  $\gamma \times B$ , known as the Larmor frequency
- Spin Alignment Occurs at the Rate, T1

Center for Cognitive Neurosciences

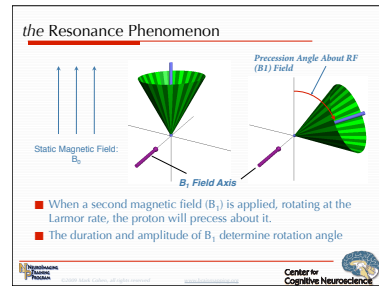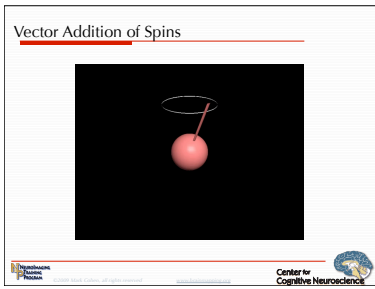


If an additional magnetic field is applied the nuclei will precess about their vector sum. In the special case that the second field,  $B_1$ , is made to rotate at the same (Larmor) rate as the proton nuclei, the  $B_1$  field will seem stationary with respect to any individual proton. They will thus precess about  $B_1$  in a simple manner. This view of the interaction of  $B_1$  and  $B_0$  (the static field) is called the rotating frame. In the rotating frame there is no apparent precession about  $B_0$ .



Because the Larmor rate at reasonable magnetic fields is typically in the tens to hundreds of MegaHertz, and because  $B_1$  is applied just long enough to produce the desired precession away from the longitudinal axis, the  $B_1$  field is typically called a Radio Frequency, or RF pulse.

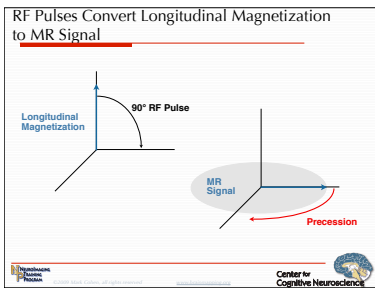


In simple terms, then, the effect of a  $90^\circ$  RF or  $B_1$  pulse is to convert any longitudinal magnetization into detectable signal.

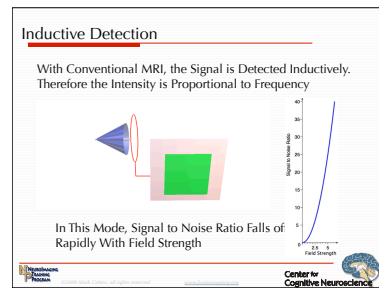


The rotating magnetic field from an ensemble of protons is usually detected by "induction". The time varying magnetic field from the spinning protons results, by Lenz's law, in the creation of a time-varying electric field. This will result in the flow of electrical current in any conductor, such as an antenna. The strength of the electrical field is proportional to the Rate of change of the magnetic field. Thus, the detection sensitivity goes up as the magnetic field is increased. For this reason, the signal to noise ratio is a quadratic function of Magnetic field strength, making inductive pickup of MRI signals at very low field effectively impossible.

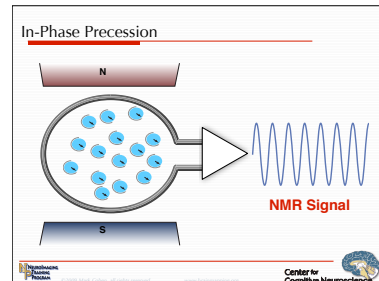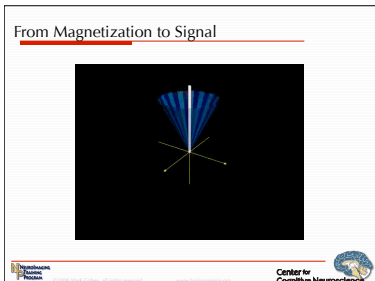


After an RF pulse, the phases of the individual spins will instantaneously be the same - that is, the spins will be precessing "in-phase" and their transverse components will add. This gives rise to a rotating magnetic field that can be detected easily because a time-varying magnetic field creates an equivalent electrical field. If a conductor (antenna) is placed in the vicinity, the electrical field will create a current that can be amplified. Generally, the signal will be a sinusoid at the Larmor frequency.

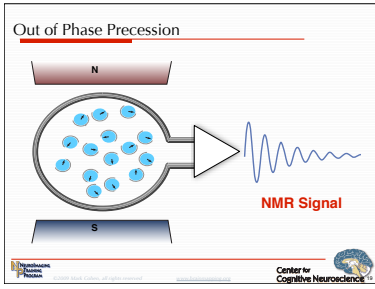




The MR signal, once formed, decays quickly. This is because the spins will precess in phase only to the extent that they precess at precisely the same frequency. Even tiny variations in field cause the dephasing to take place over a couple of tenths of a second. These field variations arise from many factors. For example, we cannot manufacture a magnet field that is homogeneous to more than a few parts per million. Further, the individual spins affect one another, causing local inhomogeneities, and different tissues magnetize to different degrees.

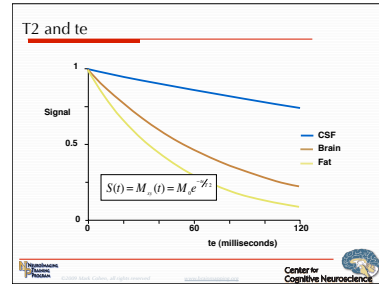


The MR signal decays exponentially as a result of spin dephasing. The decay is governed by a time constant, T<sub>2</sub>, that is tissue specific. In the body, the T<sub>2</sub> ranges from a few hundred microseconds in bone to hundreds of milliseconds in water (CSF, urine). The adjustable parameter, te, specifies the time after signal excitation at which the data are collected. Longer te's result in less signal overall, but also in increased contrast between tissues with differing T<sub>2</sub>.


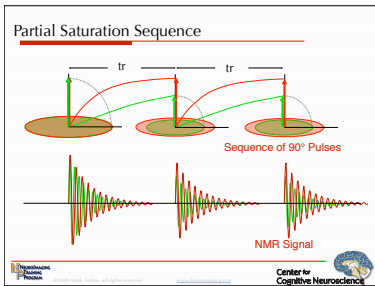


Repeating the excitation and data collection at a fixed interval, tr, introduces contrast based on T<sub>1</sub> differences. This is because if the tr is short compared to T<sub>1</sub>, the spins do not regain their equilibrium state. Each time that a 90° pulse is applied, all of the longitudinal magnetization is lost, and starts again from zero. If a tissue has a short T<sub>1</sub>, more magnetization is regained, and more signal is created with each 90° pulse. The signal from longer T<sub>1</sub> tissues is weaker, however.

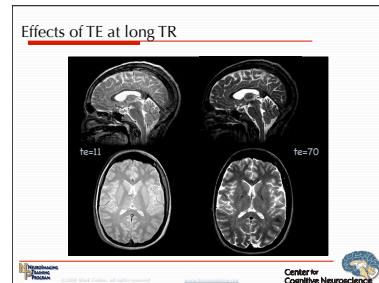
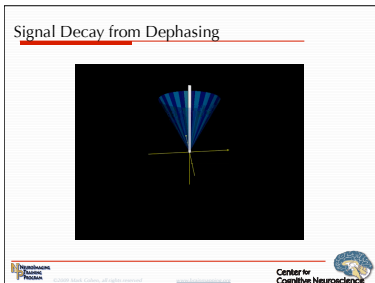
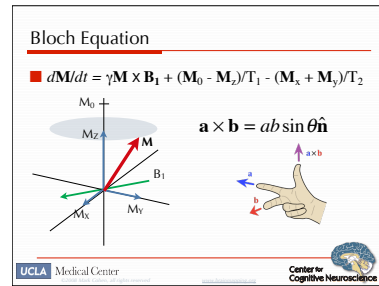


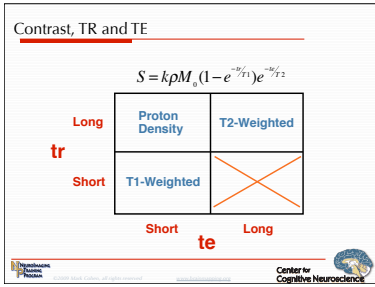




This chart suggests the four principle variations in contrast used in MR imaging. If we consider long tr to be  $tr \gg T_1$ , the images will have little T1 contrast, as all tissues will magnetize almost fully. If  $tr \ll T_1$  the images will have T1 contrast. Likewise,  $te \ll T_2$  minimizes the contrast differences from T2 and  $te \gg T_2$  maximizes T2 contrast. The equation shows this in analytic form. When tr is long, and te is short, the major contrast determinant would be the density of protons per unit volume. Images with mixed T1 and T2 contrast are typically poor, as shown in the next slide.

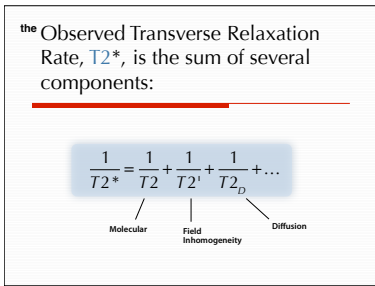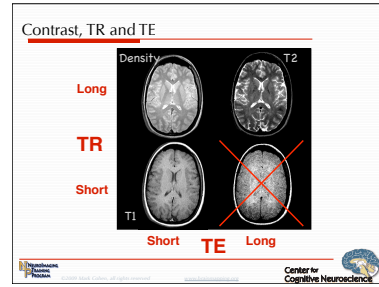

Many physical effects result in dephasing. The process that dephases the spins most quickly dominates the contrast in the final images, thus the net observed T2 ( $T_2^*$ , "T2-star") is shorter than any of the individual T2's. Most MR imaging is concerned chiefly with three dephasing effects: molecular interactions (T2), local variations in field from, for example, tissue boundaries (T2') and dephasing from molecular motion from diffusion (T2D).

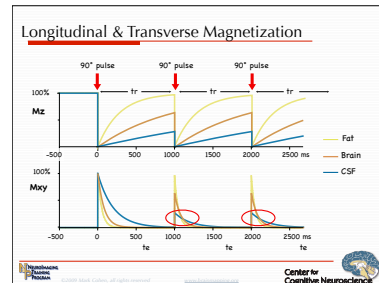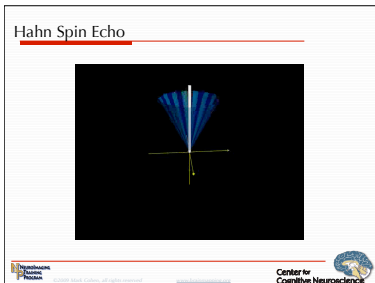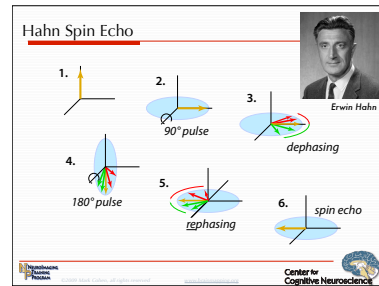

This graph shows the combined effects of T1 and T2 on the longitudinal and transverse magnetization during repeated excitation pulses at fixed tr. This model was made with a tr of 1 second, and with typical T1 and T2 of head tissues. Notice that the signals for brain, CSF and fat all cross over with a te of about 70 msec at this tr. The result is that there is no contrast between these tissues. If te is made very short the T1 effects dominate. If te is made very long, the T2 effects dominate, but the signal will be VERY weak.

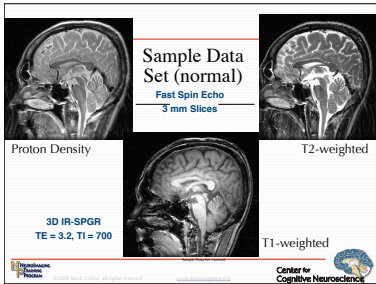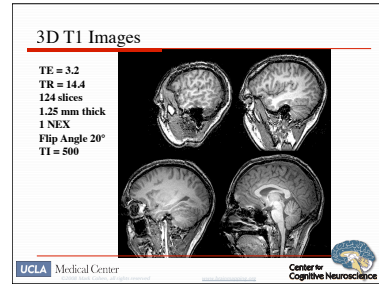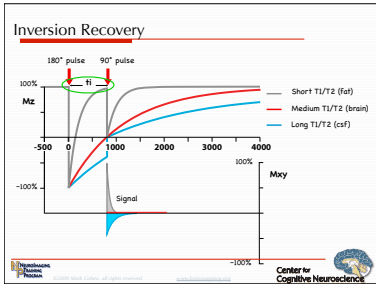

### MR Formulæ

Spin Echo Signal =  $k\rho M_0(1 - e^{-tr/T1})e^{-te/T2}$


Inversion Recovery Signal =  $k\rho M_0(1 - 2e^{-ti/T1} + e^{-tr/T1})e^{-te/T2}$

$\rho$  is the proton density  
 $k$  represents instrument effects

The "Bloch" Equation:

$$dM/dt = \gamma M \times B_1 + (M_0 - M_x)/T_1 - (M_y + jM_z)/T_2$$

Center for Cognitive Neurosciences

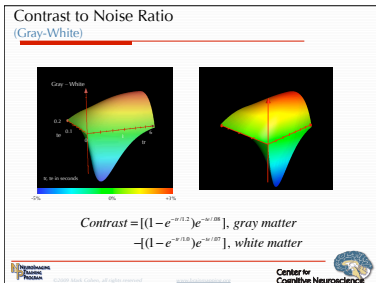


- ### MRI Contrast Summary
- Pulses of Rotating Magnetic Fields (RF) Convert Nuclear Magnetization to Signal
  - RF Pulses Add Energy by Displacing Longitudinal Equilibrium
  - Contribution of Intrinsic Tissue Properties T1 and T2 Manipulated by Experimenter controlled timing: tr and te respectively.
  - Typically,  $0.05 < T1 < 3s$  and  $0.005 < T2 < .3s$  for body tissues.
- ...next: Making and Image
- Center for Cognitive Neurosciences

The graphs above plot the difference in signal (the "contrast") between grey and white matter as a function of tr and te. Notice that there is a broad range of combined tr and te where the signals will be isointense (no contrast). The blue regions of the graphs are regions of high T1 contrast, and the red regions have high T2 contrast.


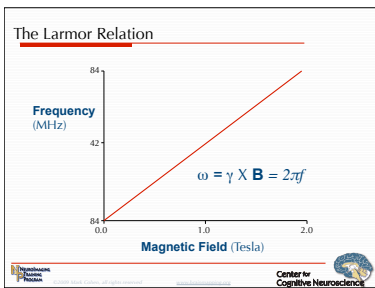
### The Plan

- The Magnetic Resonance Phenomenon & Contrast
- Spatial Encoding
- The "Pulse Sequence" Rules Everything
- Functional MRI
- Image Quality and Artifacts




### MR Spatial Encoding

- Most MR Spatial encoding is based on a single concept:
  - If the **Magnetic Field** varies by location, the **MR Frequency** will vary by location,
  - Therefore: We can *find the location by measuring the frequency.*
- Newer methods (e.g., GRAPPA, SENSE) also find location by signal strength as a function of antenna distance





Given what we know about signal intensity and contrast, our next step is to make an image. This implies finding the location of the signal in 3 dimensional space. To do this we rely on the Larmor relation that the precession frequency is proportional to magnetic field. The trick will be to cause the magnetic field to vary with position.

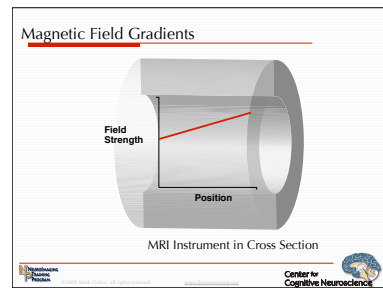

We create spatially varying ("gradient") magnetic fields within the bore of the imaging magnet through the use of physical "gradient" coils. The spatial variations are typically small compared to the strong imaging magnetic field. These cause the frequency,


### Frequency Selective Excitation

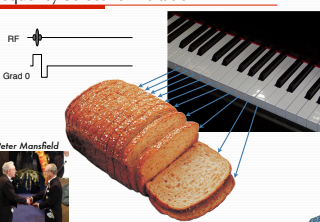


RF

Grad 0

Sir Peter Mansfield



Slice selective excitation is an invention of Nobel-laureate, Sir Peter Mansfield. He noted that if an excitation RF pulse was applied at a fixed frequency in the presence of a gradient field, only spins at the corresponding location in the field would be affected by the pulse and become rotated into the transverse plane (excited). The position corresponds to frequency and can be selected easily.

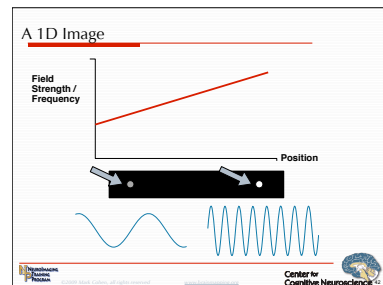


The frequency of spin precession depends on the strength of the local magnetic field. If the field varies with position, so will the frequency. In this diagram, a stronger NMR signal appears in a location having a higher magnetic field.

### Interlude: The Fourier Transform

- Any (continuously differentiable) function may be represented as the sum of a set of sinusoids of varied amplitude and frequency.
- The **Fourier Transform (FT)** is "complex". In general it produces Real (cosine) and Imaginary (sine) components.
- The Fourier transform is computed easily using the Fast Fourier Transform:
  - MATLAB: `>> F = fft(x)`

### Fourier Sum example

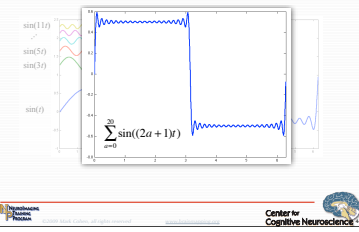

At left are 6 sinusoids at frequencies of 1, 3, 5, 7, 9 and 11 cycles, weighted by a factor determined by the Fourier transform. At right is the sum of these sinusoids. In the bottom figure, the first 20 odd frequencies are added to make a near perfect 'square' function

### The Fourier Transform

- The Fourier transform  $\mathcal{F}(\omega)$ , of a function  $f(t)$ , is formed from the products of  $f(t)$  with sinusoids at all frequencies.

The Fourier transform is a function of frequency,  $\omega$


Yikes!  $i = \sqrt{-1}$

$$\mathcal{F}(\omega) = \int_{-\infty}^{\infty} f(t) [\cos(2\pi\omega t) - i \sin(2\pi\omega t)] dt$$


$$= \int_{-\infty}^{\infty} f(t) e^{-2\pi i \omega t} dt$$

This amazing result is from Euler  $e^{ix} = \cos(x) + i \sin(x)$


### The Fourier Transform is: the Product of a function and a Series of Sines and Cosines

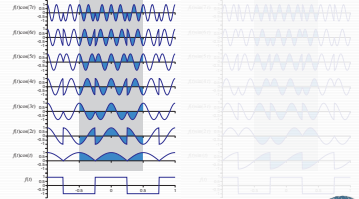


### Typical Representation of Fourier Transform

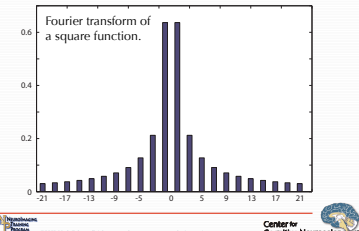


### Objects with Finite Width

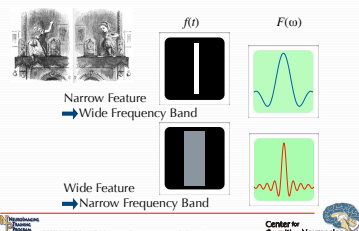




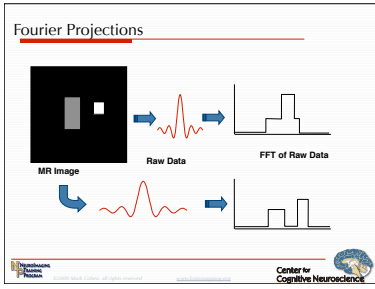


After a slice has been excited, we can apply a field gradient in a different direction, causing the spin frequencies to vary along this second axis. The Fourier transform of the signal represents calculates the intensity of the signal at each frequency. As frequency depends on position, we can interpret this as the intensity as a function of location along the same axis as the gradient field was applied. We can readily apply gradient fields along any axis, and thus collect these "projections" in any orientation.

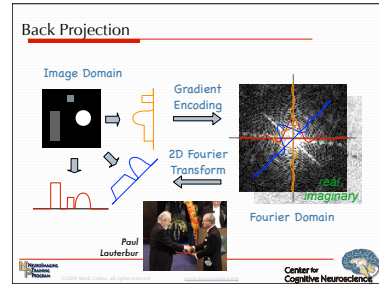


Here we show the projected and transformed signal along three different axes. We can project the 'untransformed' data lines into a raw data space according to the direction of the gradient field. The combination of all of these radial projections can look a bit like rings of waves in a puddle. The 2 dimensional Fourier transform of this combined raw data is the image. A 2D Fourier transform amounts simply to calculating the Fourier transform of the data set line by line, then calculating the Fourier transform of each column in this intermediate data set.

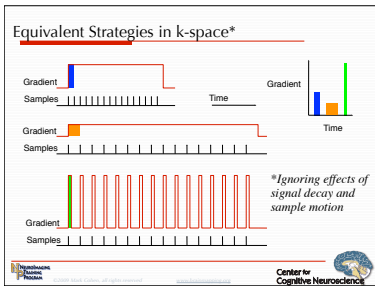


Crucially, if we ignore all of the contrast effects considered above, the only thing that causes the intensity of the MR signal to change is the effect of the applied magnetic field gradients. While the signal evolves in a continuous manner, in practise, it is sampled discretely and digitized. In the figure, the circles on stems represent sampling time points. In the time between samples (incidentally, this is called the "dwell time") the differences in precession frequency of spins in different locations causes a phase difference to accumulate as a function of the product of the frequency difference between these locations and the dwell time. Thus, leaving the gradients on for a longer time at a shorter amplitude. From the perspective of the difference between each sample time point there is no difference. In fact, we don't have to leave the gradients on for the duration of the dwell interval and instead, we may pulse the gradients at very high amplitude and get the same effect.

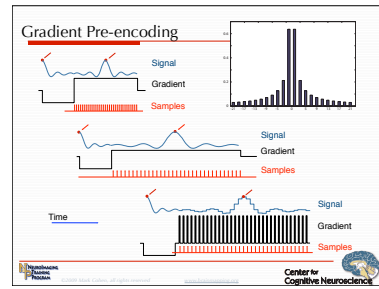

The pulsed gradient method (bottom) causes the signal to make discrete jumps in intensity. If the gradients are left on continuously, but the signal is sampled discretely, we collect identical data. Note in these figures that the spins are typically "pre-encoded" by applying the field gradients in their opposite sense, before collected the signal. This ensures that all of the spins will be in-phase at the center of the readout period.

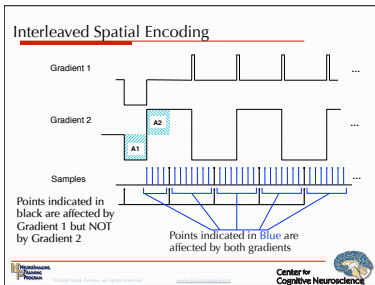

As a gradient is applied for a time period (A1) the spins in different location go out of phase with respect to one another. When the gradient is applied for the same duration in the opposite sense - reversing which side of the instrument is at higher field and which is lower - the spins go back into phase. Considering only Gradient 2 above, at the time points indicated in black, the spins are always in phase. Alternating the sign of Gradient 2 makes it invisible at these time points. However, if an orthogonal gradient, Gradient 1 is applied between these oscillations, it will cause a phase accumulation along the Gradient 1 direction. For example, if Gradient 1 is directed along the vertical (Y) axis of the magnet, the time points indicated in black represent a line with spatial encoding along this axis only. The Fourier transform of this, of course, is a signal intensity projection along the Y axis. The points indicated in blue are affected by both gradients and fill in the other points in the 2D raw data space. In this way, we can interleave the spatial encoding in the X and Y axes creating a 2D image all at once. This scheme is known as Echo Planar Imaging or EPI and is the fastest practical method to form MR images at this time.

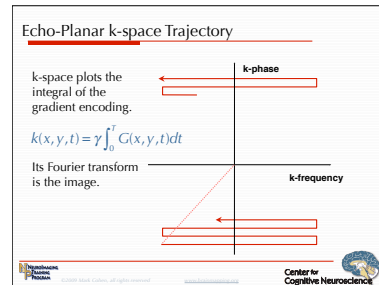

If we examine the locations of the raw data points that are collected in EPI we see that they are collected line-by-line in alternating directions. The brief pulses of Gradient 1 move the locations from one line to the next. The net effect is to make a raster pattern in the raw data space. This raw data space is known as "k-space" and is the 2D Fourier transform of the image. As long as there are no other things changing the signal, the order in which the data are placed in k-space makes no difference. Since the amount of spatial encoding between spins is determined by the product of the gradient amplitude and duration we generalize this as the integral of the gradient-time product. This product becomes the phase difference between spin locations; k-space is sometimes referred to as phase space for this reason.

### Properties of K-Space

- Increasing K values Represent Higher Spatial Frequencies, thus Higher Resolution
- Finer Grain Sampling Results in Wider FOV
- Reflections Across  $K=0$  are approximate Complex Conjugates

Center for Cognitive Neuroscience

### Raw Data Symmetry

Center for Cognitive Neuroscience

### k-space conjugate symmetry

For a Stationary Object, in a Homogeneous Field:  $S(kx, ky) = S(kx, -ky)$

where  $S(kx, ky)$  is the signal at  $(kx, ky)$ .

Example: if  $S(kx, ky) = a + ib$ , then  $S(kx, -ky) = (a + ib) = a - ib$ .

Center for Cognitive Neuroscience

### Conventional Spatial Encoding

Center for Cognitive Neuroscience

Before echo-planar imaging became possible, a different encoding scheme was used that collected only one line of k-space with each  $tr$ . The position along the  $y$ -axis of k-space is established by applying a brief pulse of the  $Y$ -gradient (Grad 1) before collecting the data in the presence of the  $X$ -gradient. With each  $tr$ , the area of the  $Y$ -gradient pulse is increased. The diagram above also shows the selective excitation pulse, wherein an RF pulse is transmitted in the presence of a gradient directed along the  $z$ -axis (Grad 0). Patterns of this kind that indicate the timing of events in MRI are known as pulse sequences.

### Conventional K-Space Trajectory

Center for Cognitive Neuroscience

The traditional line-by-line manner of forming MRI images results in the k-space filling pattern shown above. Note that time does not appear explicitly in the k-space diagram. In the conventional imaging case, the total imaging time is determined by the number of lines of resolution in the  $Y$  or phase encoding axis and by the time between lines, which is  $tr$ .


### Spiral

Center for Cognitive Neuroscience

The order in which k-space is populated is somewhat arbitrary. One filling pattern, "Spiral Scanning" has significant popularity. The idea is to make a spiral trajectory in k-space. It is easy to see that this requires a gradient pattern in  $X$  and  $Y$  that is a slowing sinusoid.


### The Plan

- The Magnetic Resonance Phenomenon & Contrast
- Spatial Encoding
- The "Pulse Sequence" Rules Everything
- Functional MRI
- Image Quality and Artifacts




### A Pulse Sequence Controls

- Slice Location
- Slice Orientation
- Slice Thickness
- Number of Slices
- Resolution (FOV and Matrix)
- Contrast (TR, TE, TI, Flip Angle, Diffusion, etc...)
- Artifact Correction (Saturation Pulses, Flow Comp, Fat Suppression, etc...)



The pulse sequence ultimately controls a host of factors. The user seldom interacts directly with the gradient and RF pulse waveform controls, but instead enters parameters for acquisition such as the slice locations resolution features and contrast controls.


### Crib Sheet

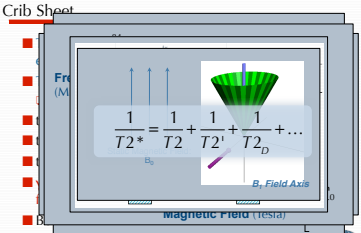



$\frac{1}{T2^*} = \frac{1}{T2} + \frac{1}{T2'} + \frac{1}{T2D} + \dots$

magnetic field (Tesla)

B1 - Rotating Magnetic Field





### Conventional Spatial Encoding

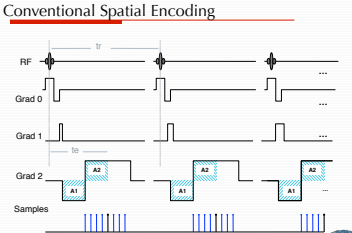Before echo-planar imaging became possible, a different encoding scheme was used that collected only one line of k-space with each tr. The position along the y-axis of k-space is established by applying a brief pulse of the Y-gradient (Grad 1) before collecting the data in the presence of the X-gradient. With each tr, the area of the Y-gradient pulse is increased. The diagram above also shows the selective excitation pulse, wherein an RF pulse is transmitted in the presence of a gradient directed along the z-axis (Grad 0). Patterns of this kind that indicate the timing of events in MRI are known as pulse sequences.


### Contrast Encoding

Before echo-planar imaging became possible, a different encoding scheme was used that collected only one line of k-space with each tr. The position along the y-axis of k-space is established by applying a brief pulse of the Y-gradient (Grad 1) before collecting the data in the presence of the X-gradient. With each tr, the area of the Y-gradient pulse is increased. The diagram above also shows the selective excitation pulse, wherein an RF pulse is transmitted in the presence of a gradient directed along the z-axis (Grad 0). Patterns of this kind that indicate the timing of events in MRI are known as pulse sequences.
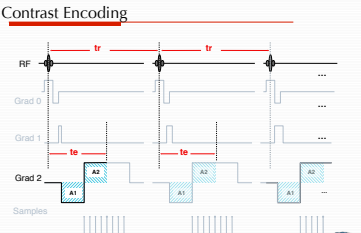
### The Plan

- The Magnetic Resonance Phenomenon
- MR Signal Formation and Contrast
- Spatial Encoding
- The "Pulse Sequence"
- Reduced Flip Angle & Steady State Approaches
- Diffusion



## Reduced Flip Angle Imaging

### Outline

- Determinants of Imaging Time
- TR, Saturation and Image Quality
- Reduced Flip Angle Techniques
  - FLASH (=SPGR)
  - FISP (=GRASS)
- Gradient Echoes
- Applications of Shallow Flip Imaging
- Ultra-Fast Imaging

## Determinants of Imaging Time

Scan Time =

$$\text{Repetition Time (TR)} \times \text{Number of Phase Encodes} \times \text{NEX (Averages)} \times \text{Number of 3D Steps}$$


## TR and Image Quality

### Reduced TR Yields:


- Decreased Scan Time
- Increased T1 Contrast
- Reduced (Useable) T2 Contrast
- Reduced Signal to Noise Ratio
- Increased Power Deposition
- Reduced Slice Coverage

## Signal and Flip Angle

Small Flip Angle

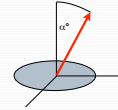


Large Flip Angle

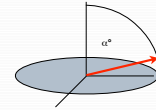


## Small and Large Flip Angle

Loss of Longitudinal Magnetization

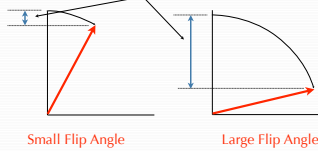

## Flip Angle and TR/T1

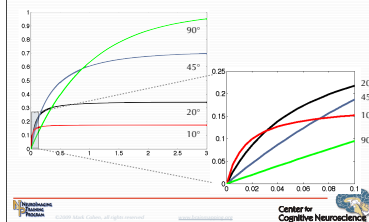


## Contrast and Flip Angle


Large Flip Angles	Short	Long
	Long	Proton Density
Short	T1 Weighted	

Small Flip Angles	Short	Long
	Long	Proton Density
Short	Proton Density	T2* Weighted

## Spin Echo Sequence

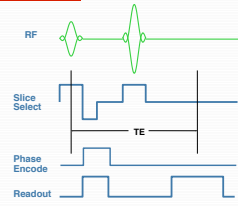

## A 180° Pulse is not used in FLASH imaging

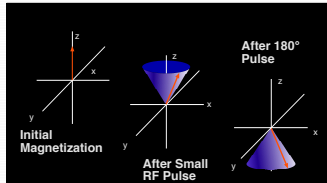

## T2 and T2\*


- T2: Transverse Magnetization Decay from Spin-Spin Interactions
- T2\*: Transverse Magnetization Decay from Local Magnetic Field Variations





## Magnetic Susceptibility

The Extent to Which a Substance Becomes "MAGNETIZED" when Placed Within a Magnetic Field



## Magnetic Susceptibility

Objects with Susceptibility Different than Air Distort the Magnetic Field

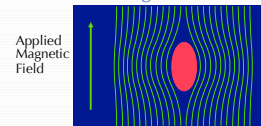

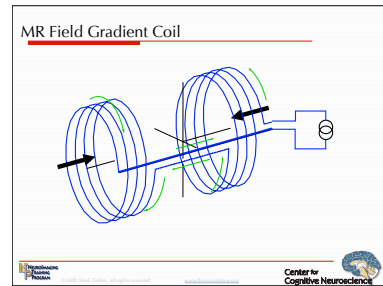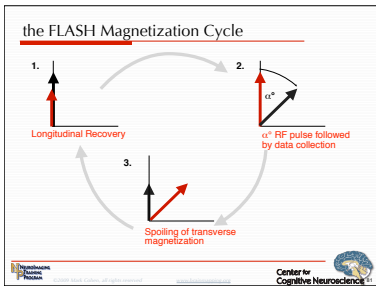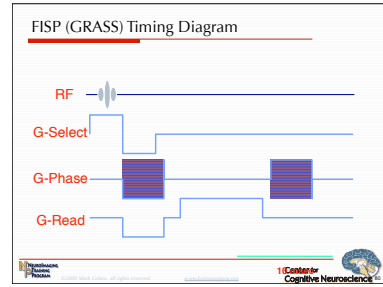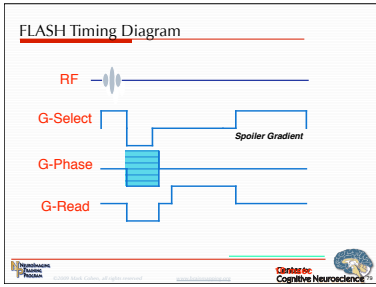

Creating the field gradients requires passing electrical current through large coils (Gradient Coils). In MR imaging systems several hundred Amperes of current are required to create large enough gradients. The pulse sequences generally require that the gradients are turned on and off as rapidly as possible. Doing so requires very high Voltages because of the coil inductance.

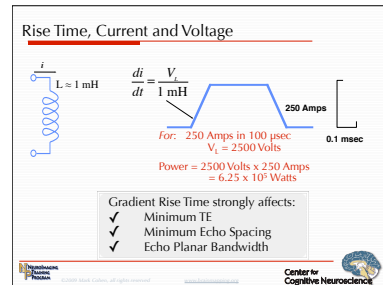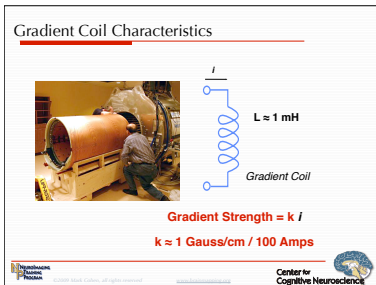


### An "Equation" in Resolution

Because MR is an emission modality the temporal resolution, spatial resolution and contrast are inter-dependent:

$$\text{Signal} = kB_0(\text{voxel size})\sqrt{\text{imaging time}} - \text{contrast}$$


where  $B_0$  is the field strength.

### Why is MRI So Noisy?

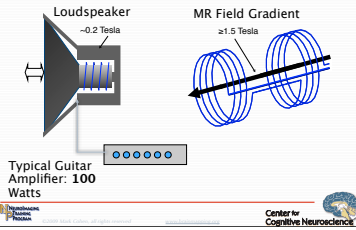

The construction of a gradient coil is very similar to the construction of a loudspeaker. Both utilize a coil of wire to produce a time-varying magnetic field in the presence of a large stationary field. The magnetic field in the coil experiences forces against the stationary field that tend to make it move. However, the magnitudes of all of the forces are much greater in MRI. Despite great efforts to reduce it, MRI is necessarily extremely noisy.

### Why is MRI So Noisy?

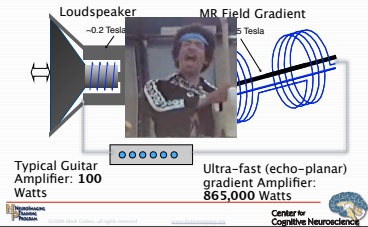

The construction of a gradient coil is very similar to the construction of a loudspeaker. Both utilize a coil of wire to produce a time-varying magnetic field in the presence of a large stationary field. The magnetic field in the coil experiences forces against the stationary field that tend to make it move. However, the magnitudes of all of the forces are much greater in MRI. Despite great efforts to reduce it, MRI is necessarily extremely noisy.

### Spatial Encoding Summary

- Spatial Encoding and Contrast are Linked through the Pulse Sequence;
- The MRI Raw Data are the 2D Fourier Transform of the Final Image (usually the Magnitude Transform);
- Spatial Encoding is Added through Gradient Coils;
- Hermite Symmetry of the Raw Data may be used to Reduce Scan Times;
- Literally Hundreds of Pulse Sequences are in Common Use,

### The Plan

- The Magnetic Resonance Phenomenon & Contrast
- Spatial Encoding
- The "Pulse Sequence" Rules Everything
- Functional MRI
- Image Quality and Artifacts


### Brain "Activation" Leads to:

CBF	Increased	+ΔR1
CBV	Increased	+ΔR2 (C+)
O <sub>2</sub> Utilization	Increased slightly?	
Venous [O <sub>2</sub> ]	Increased	-ΔR2* ← "BOLD"
Glucose Utilization	Increased	? Lactate

R1=1/T1  
R2=1/T2


There are many biophysical effects associated with increased neuronal electrical activity. At least five coupling events are known to be observable during fMRI, though most studies at present rely on detecting the change in venous oxygen concentration which follows increases in synaptic activity.

### Signal Losses from Spin Dephasing

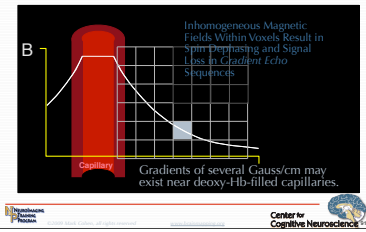

MRI signal depends on the coherence of individual proton spins. In the presence of magnetic field distortions, the MR signal is reduced. Deoxygenated blood is paramagnetic and results in such field distortions in the immediate vicinity of the capillaries.


### Why Does Venous O<sub>2</sub> Increase? ..

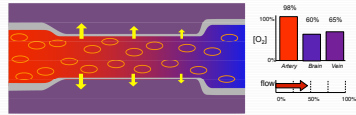

Under normal conditions oxygen diffuses down its concentration gradient from the capillary to the brain parenchyma


This series of slides outlines the dominant hypothesis as to why the BOLD signal is increased when brain metabolic demand is increased. Some believe, however, that this demand is driven not by oxygen but by glucose. The effects would be the same, however.

### Why Does Venous O<sub>2</sub> Increase? ..

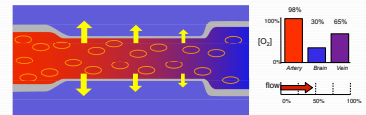


As the brain becomes more active, the oxygen consumption increases, increasing the transmural oxygen gradient.

### Why Does Venous O<sub>2</sub> Increase? ..

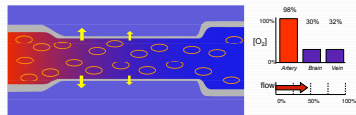


As oxygen flows across the capillary lumen it is depleted in the capillary and no further oxygen can be delivered

### Why Does Venous O<sub>2</sub> Increase? ..

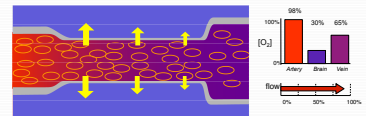

The vascular system responds by increasing blood flow so that more oxygenated blood is available throughout the capillary

### Why Does Venous O<sub>2</sub> Increase? ..

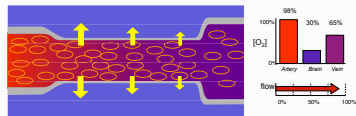

Because the blood flow is increased more oxygenated blood passes into the venous end of the capillary



### Why Does Venous O<sub>2</sub> Increase? ..

Because the blood flow is increased more oxygenated blood passes into the venous end of the capillary

Center for Cognitive Neuroscience

### BOLD Contrast & Field Strength

- BOLD Contrast arises from susceptibility differences
- The absolute field distortion (from BOLD) is proportional to the magnetic field strength
- The absolute change in MRI signal is proportional to *both* the field distortion and the signal strength.

BOLD *should* go as  $kB_0^2$

Center for Cognitive Neuroscience

### fMRI

explores intensity variations in MR signal

intensity variations reflect venous [O<sub>2</sub>]

Center for Cognitive Neuroscience

### Gradient-Recalled Echo

Photic Stimulation – GE Images

Center for Cognitive Neuroscience

fMRI Signal increases relatively slowly (several seconds), and returns to baseline slowly. No exogenous contrast needed

### Ken Kwong

Inversion Recovery  
TE=42 IR=3000  
TI = 1100  
Thickness=10

Center for Cognitive Neuroscience


A second method is sensitive to blood flow changes, and demonstrates similar results.

### Inversion Recovery


Photic Stimulation – IR Images

Center for Cognitive Neuroscience


Note that the flow signal continues to increase for a longer time. This is due to complex effects of perfusion and proton exchange, as well as to possible vascular dynamic effects

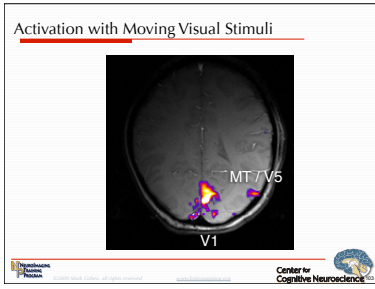


One of the more interesting and important features of fMRI is its ability to collect multiple experiments on the same subject to tease apart the responsiveness of various brain regions to complex stimuli. In this study, Roger Tootell looked at the fMRI signal change when a subject observed a visual stimulus containing moving vertical bars. Multiple brain areas showed signal increases, but his hypothesis was that the detection of the bars and of their motion were represented separately in the brain.

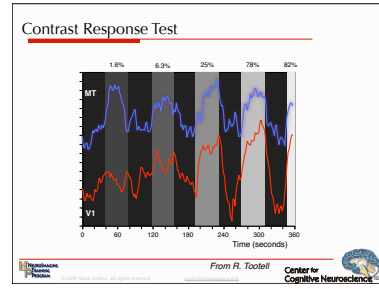


To test this hypothesis, Tootell looked at the fMRI signal as a function of contrast between the moving bars and their background. As the contrast was increased, there was a parametric increase in the signal response within area V1 (the central blob in the prior slide) whereas the more lateral extrastriate blob responded similarly regardless of contrast.

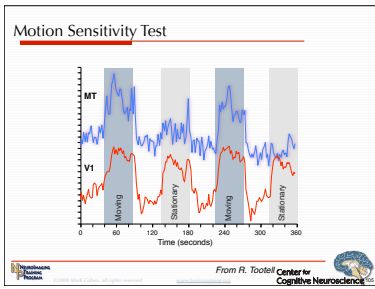


When Tootell instead observed the responses to lines that were either moving or stationary, the signal in V1 was of essentially equal amplitude whether or not the bars were moving, whereas the signal changes were much greater in the extrastriate area, MT, when the bars were moving. This is a double dissociation experiment that was able to functionally identify area MT.

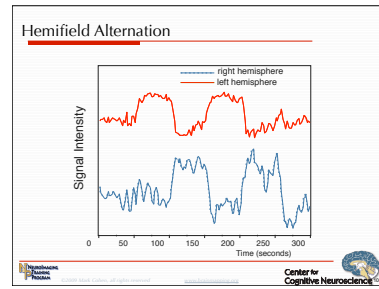


If stimuli are modulated relatively slowly, the BOLD fMRI response follows the stimulus timing rather faithfully.

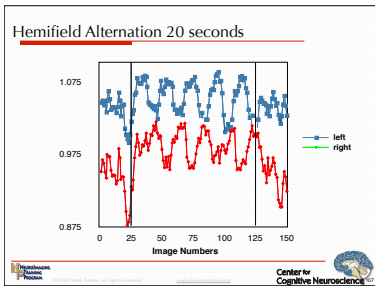


When stimuli are varied quickly, the BOLD signal changes are much reduced and do not look a great deal like the stimulus timing.

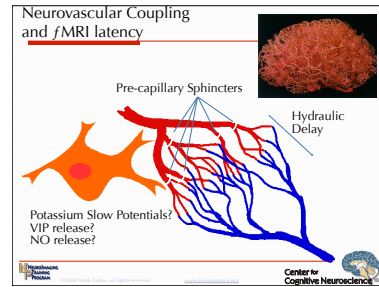


The fMRI response takes some time to occur. Much of this delay is thought to be related to the signaling from neurons to the vascular system. At this time it is not entirely clear how this signaling takes place, though there are many candidate signaling molecules and mechanisms that are probably all involved to varying degrees.

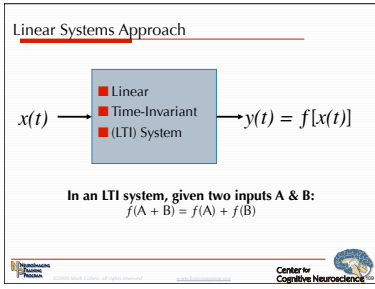


Almost all fMRI analysis is now performed under the rubric of linear systems analysis. This is because the properties of linear systems are very well-worked out and the toolsets available to analyze them are mature and relatively simple. A linear, time-invariant system is defined as one which obeys the property that the response of the system to two inputs presented together is the same as the sum of the responses to those two inputs presented individually. One of the most important properties of a linear system is that by measuring its response to a standard input - usually an instantaneous transient input, we can calculate and predict its response to any input through the use of convolution. Specifically, the predicted output, Y, to any stimulus, X is the convolution of X with the systems response, h, to an impulse:  $Y = X * h$ , where \* is used to indicate convolution.

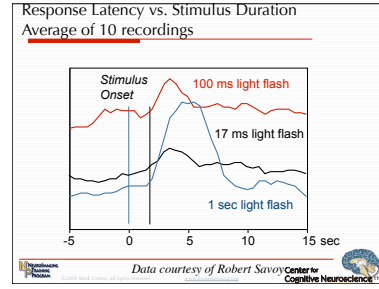



To capture the power of linear systems analysis for fMRI data it is useful to estimate the form of the brain's response to an impulse stimulus. These data of Robert Savoy show the response in V1 to brief flashes of light, which is a reasonable approximation to an impulse stimulus.

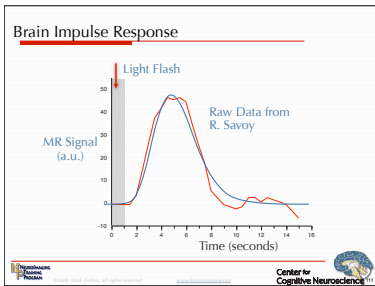


Using the observed V1 impulse response one popularly used model of the brain impulse response is the Gamma variate waveshape, shown here in blue. The Gamma function is a relatively function whose shape is very similar to the observed brain impulse data.

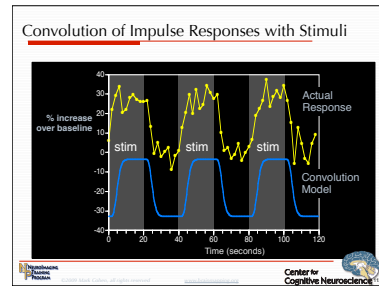


The yellow points indicate the signal intensity observed in V1 in response to an intense flashing light stimulus that was turned on and off every twenty seconds. Convoluting the time course of the stimulus with the Gamma function discussed above results in the estimate for brain activity shown in blue. Visual comparison of the two curves suggests that the convolution model is a reasonably good approximation to the observed response. The canonical means of detecting brain activation by fMRI is to model the brain response to the stimulus and to search for brain regions that behave similarly to the model. Thus, the convolution approach based on linear systems theory is a powerful tool for this work.

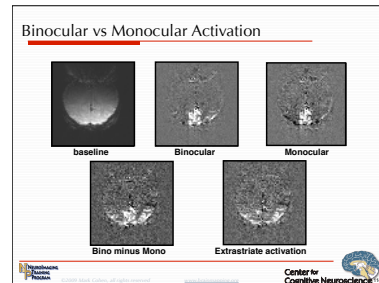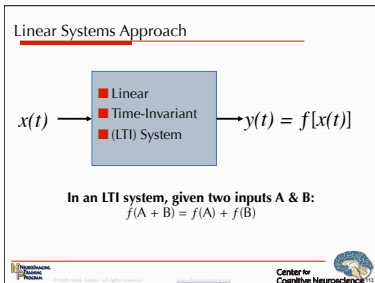

Given the importance of the linear systems analysis, it is instructive to test whether the brain fMRI signal response is linear at all. These data from Ken Kwong show the brain responses to visual stimulation to both eyes simultaneously and to one eye alone. Because the primary visual cortex is known to be independent for each eye, we would expect the difference between the binocular and monocular stimulation to be the same as the response to the other eye alone. The data however show that this is simply not the case. Thus our assumption of linearity is strongly violated.


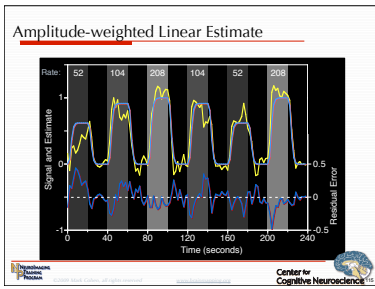


Various degrees of non-linearity can be incorporated into the linear systems analysis in a principled manner. Here, the estimated BOLD response was modeled as increasing in amplitude with the log of the rate at which subjects performed a finger motion task. Notice that the response estimates closely fit the data.  
 People using fMRI for their data collections should be cognizant of the assumptions made in the analysis and thoughtful of how they relate to the underlying physiology of the brain. In many cases, the assumptions of linearity will produce grossly distorted activation maps.

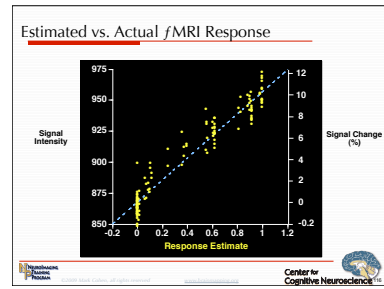

These are the data from the prior slide plotted instead as the observed signal intensity as against the modeled intensity and show that the log transform model closely fits the actual data.

- ### The Plan
- The Magnetic Resonance Phenomenon & Contrast
  - Spatial Encoding
  - The "Pulse Sequence" Rules Everything
  - Image Quality and Artifacts
  - anti-aliasing - *apologia*
  - Diffusion

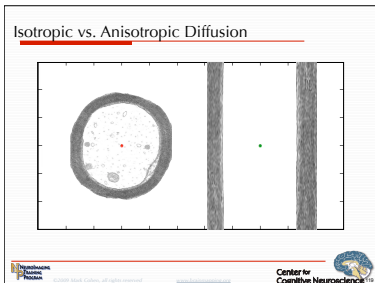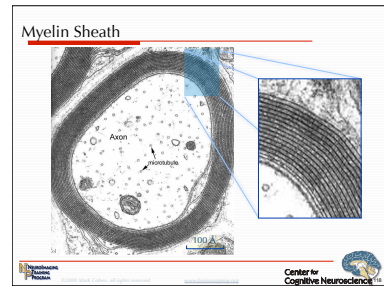

### Motion Effects on Phase






- Relative Phase is the product  $\gamma B t$  (cycles/sec/Tesla \* Tesla \* sec)
- Motion causes additional phase shifts:

$$x(t) = x_0 + vt + at^2/2 + \dots$$







$$\phi = \gamma \int_0^T G(t)x(t)dt$$

$$= \gamma \int_0^T G(t)[x_0 + vt + at^2/2 + \dots]dt$$

$$= \gamma \left[ \frac{vT^2}{4} + \frac{23at^3}{24} + \dots \right]$$

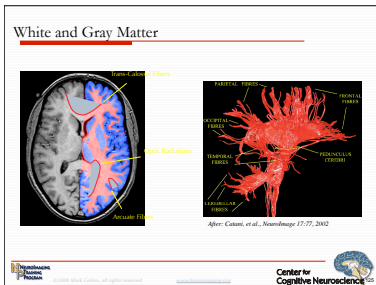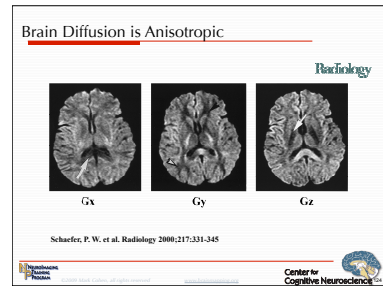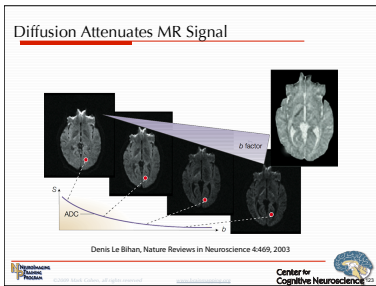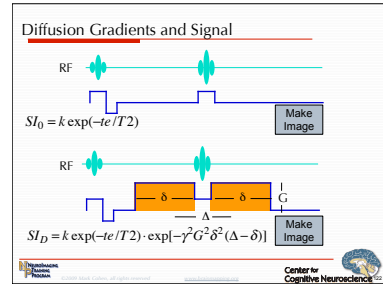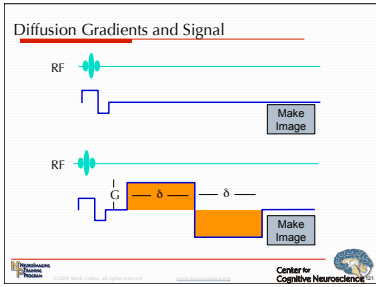


- ### The Plan
- The Magnetic Resonance Phenomenon & Contrast
  - Spatial Encoding
  - The "Pulse Sequence" Rules Everything
  - Functional MRI
  - Diffusion
  - Image Quality and Artifacts
- Center for Cognitive Neuroscience

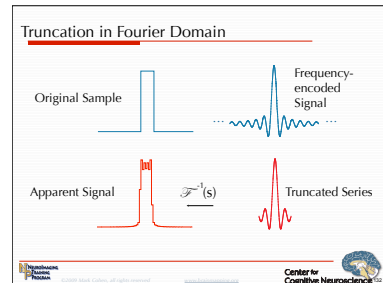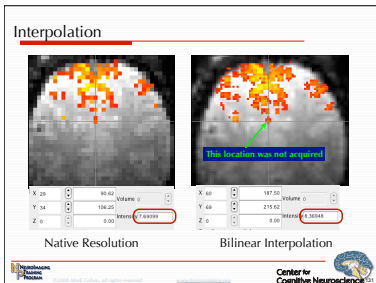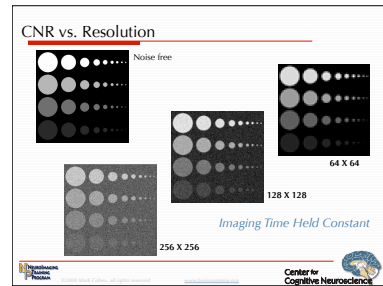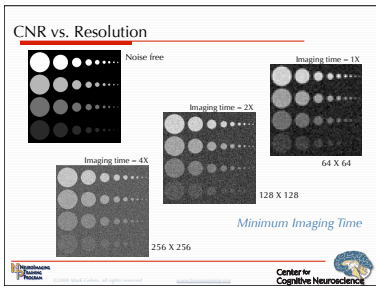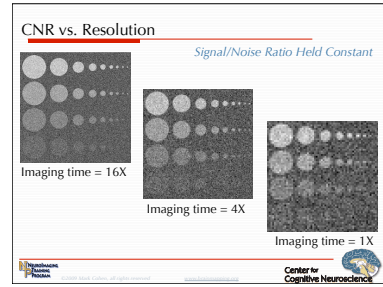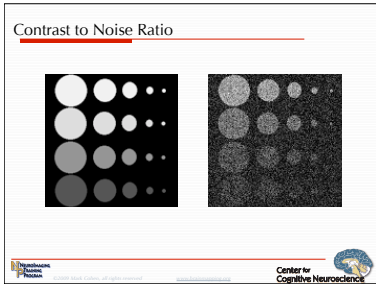


Remembering that the MRI raw data encode position by frequency, it is useful to consider factors that control the fidelity of that encoding. The Fourier transform represents the amplitudes (and phases) of the frequencies in a signal or waveform. Consider a simple rectangular signal as shown at top left. Its representation in frequency is shown at top right. Notably, the transform has infinite tails because the waveform contains energy at all frequencies.

If the waveform at upper right were the MR signal for a rectangular object, we would necessarily be able to acquire only part of it as doing otherwise would require infinite time. When we then take the Fourier transform of the truncated signal, we find a distorted representation of the original waveform. Notably, the signal now contains ringing at the edges.

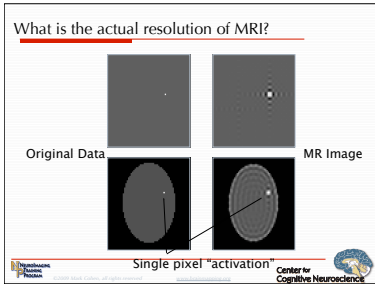

Because actual MR imaging is always based on a truncated signal collection, the resulting images are always distorted. In this simulation, we consider the images that would result from a small area of bright signal. The image representation will spread far away from the signal source as a sequence of dark and light bands. Further, a distinct, though low intensity, focus of signal will appear in the opposite X and Y locations on the image. If the bright location were a focus of brain activation that caused a BOLD signal increase, the apparent location of the activity would be represented throughout the image. An accurate statistical map based on the images would include signal decreases in the pixels next to the center of the activation.

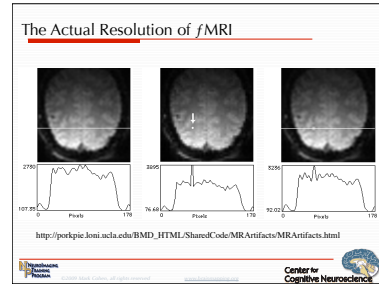

The apparent representation on an MR image of a signal point of true signal increase

### Bandwidth and Readout

- Position is encoded by FREQUENCY
- Bandwidth refers to the Frequency Difference from the center of the image to its edge:
 
$$\text{Frequency per pixel} = \frac{2 * \text{Bandwidth}}{\text{number of pixels}} = \frac{1}{\text{readout duration}}$$
- Bandwidth decreases with readout duration:
 
$$\text{Bandwidth} = \frac{\text{number of pixels}}{2 * \text{readout duration}}$$

Center for Cognitive Neuroscience

### Bandwidth and SNR

Decreasing the Bandwidth Improves SNR:  
Imaging Time is INCREASED and high frequency noise is excluded



Signal Intensity

Narrow Bandwidth

Frequency

Noise

Center for Cognitive Neuroscience

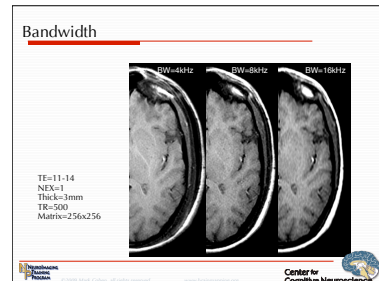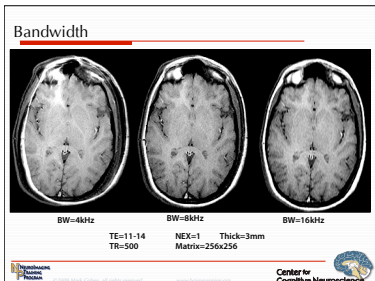

### Shape and Bandwidth

Center for Cognitive Neuroscience

### Distortions are More Severe at High Magnetic Field Strength

Variation in sample magnetization is proportional to field strength.  
High Field images lose more signal from field inhomogeneity

Center for Cognitive Neuroscience

Location is detected by the frequency of spin precession, which in turn is determined by the local strength of the magnetic field. If the field is distorted, the images will be as well. Unfortunately, the subject him or herself distorts the primary field because body tissues magnetize unequally. The images become distorted by this. Various pulse sequence factors either tend to mitigate or amplify this effect, but the combination of parameters used to create fMRI images tends to amplify these shape distortions.

### Apodization from Long Readouts

Center for Cognitive Neuroscience

If the readout period is long compared to  $T_2^*$ , the images tend to become blurred or apodized. This simulation shows the blurring that would occur with various readout durations.  $T_2^*$  in the human brain is typically around 40 msec (but depends strongly on the scanner itself) and readout periods of 50 to 80 msec are not uncommon. Thus the images are detectably blurred. This simulation considers only rectilinear k-space traversals, such as EPI. In spiral scans the blurring that results is more circular and much more complex.

### EPI Readout Durations

Center for Cognitive Neuroscience

Henceforth all of the discussion about gradient encoding and image distortion pointedly neglected T1 and T2 effects on the signal. The actual effects on the spatial localization are not insignificant, however. During the time that we are encoding the location, the signal is changing. For example, during the readout of an echo-planar or spiral scan,  $T_2^*$  effects cause the signal to decrease. In long readout sequences, such as spiral scans the signal may decrease by 80% during the actual readout.

### The Origin of Chemical Shift

Center for Cognitive Neuroscience

### Chemical Shift Artifact

If the frequency width of each pixel is less than the frequency difference between water and lipid, then water and lipid will appear in separate pixels

Center for Cognitive Neuroscience





## Chemical Shift

The Fat-Water chemical shift is about 3.5 ppm or:


Which is:

75 Hz @ 0.5 Tesla	< 1 pixel
150 Hz @ 1.0 Tesla	≈ 1 pixel
220 Hz @ 1.5 Tesla	> 1 pixel
440 Hz @ 3.0 Tesla	≈ 3.5 pixels

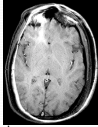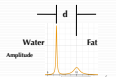


Lowering the Bandwidth/pixel increases the Chemical Shift in pixels

## Motion Artifact

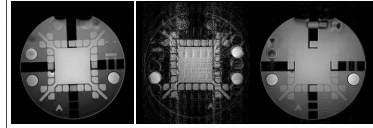


[http://artim.com/wiki/wiki.php?title=Motion\\_Correction\\_in\\_MRI](http://artim.com/wiki/wiki.php?title=Motion_Correction_in_MRI)


## Aliasing

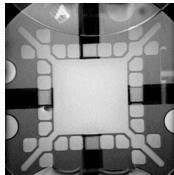


## Saturation

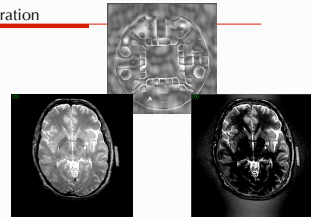

## Spikes

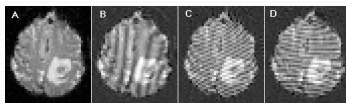

## Spikes

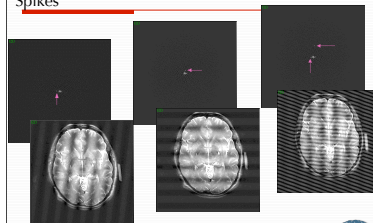

### Image Quality


- SNR is Very Limited in MRI
- Feature Detection Falls Rapidly with Loss in Contrast to Noise Ratio
- Usable Resolution is NOT the Same as Voxel Size
- Spatial Encoding Artifacts in MRI May Have Complex Appearance
- Edge Ringing and Blurring are Related to Parameters Such as Contrast

### Characterize Your Tools


- Test Statistics are Effect/Variance
- Variance includes:
  - Intra-subject (motion, attention, physiology, fatigue,...)
  - Inter-subject variance (position, morphology, performance, pathology, physiology,...)
  - Experimental Variance (uncontrolled variables, stimulation variance,...)
  - Instrument Variance
  - Site-wise Variance
  - True Random Noise

### Instrument Variation


- 1. System Instability

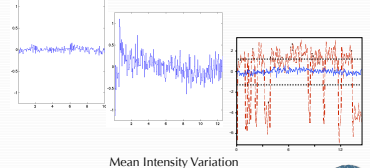

### Scanner Autocorrelation

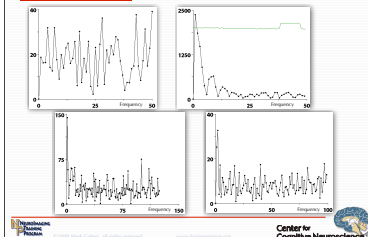


### Weisskoff Plot Weisskoff R. Magn Reson Med 36:643






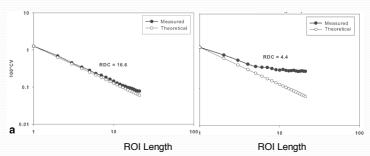


Friedman and Glover, JMRI 23:827

### Weisskoff Plot Weisskoff R. Magn Reson Med 36:643

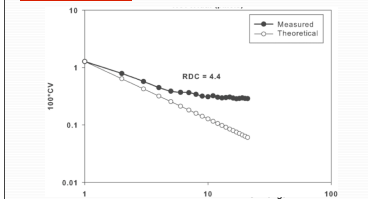

Friedman and Glover, JMRI 23:827

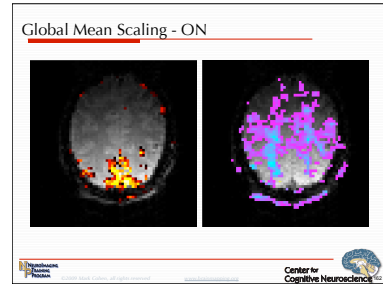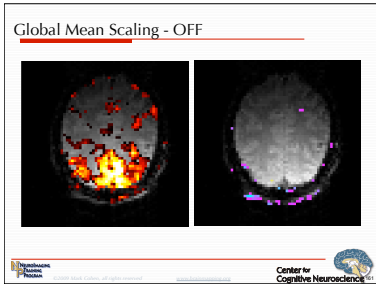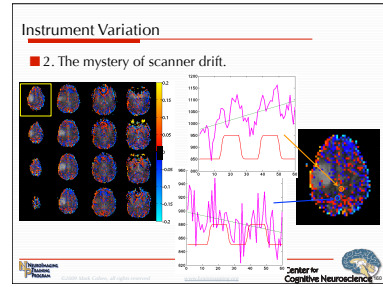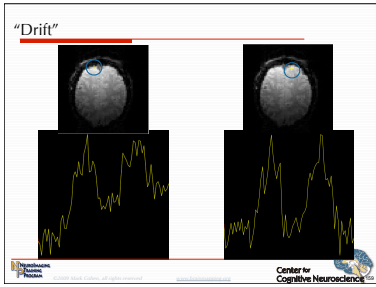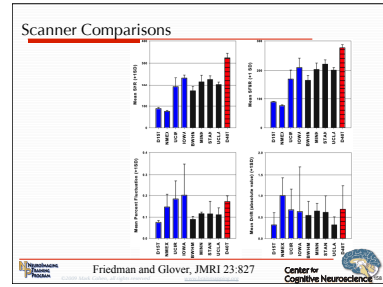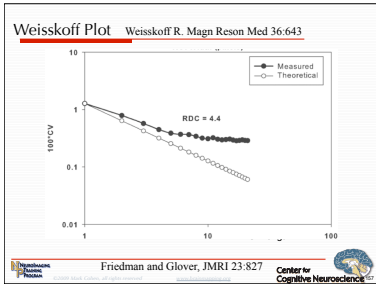


### Thermal Noise

- Noise Distribution in MRI is **Rician**
  - $Signal = \sqrt{(\Re + \sigma_s)^2 + (\Im + \sigma_s)^2}$
- Background should be **Rayleigh Noise** ( $\Re = \Im = 0$ )
  - Expected  $\mu$ :  $\sigma_s \sqrt{\frac{\pi}{2}}$  expected variance:  $\frac{4 - \pi}{2} \sigma_s^2$
  - $\left( \frac{\mu}{\sqrt{\text{variance}}} \right) = \frac{\sqrt{\pi}}{\sqrt{4 - \pi}}$
- Deviations from This Model Imply Coherent Artifacts

### Parameter Optimization



- Best BOLD contrast when  $te = T2^*$

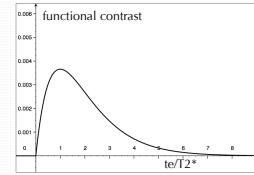

### Parameter Optimization

- Best BOLD contrast when  $te = T2^*$


$$\frac{d \ln(\text{Signal})}{d(te)} = \frac{-1}{T2^*}$$

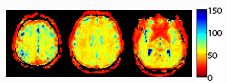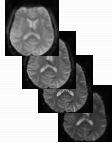

T2\* is both an *Instrument Parameter* and a *Physiological Parameter*

### Parameter Optimization - Flip Angle

- SNR per unit time is a *strong* function of  $tr/T1$

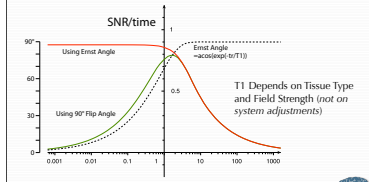

### Other Parameters...

- Voxels Should be Large enough that:
  - Thermal Noise  $\ll$  Physiological Fluctuations
- For Best Signal with Arbitrary Slice Orientation:
  - Voxels Should be Isotropic
- With Gradient Echo Scans (most BOLD):
  - Signal Falls (much) more than linearly with Slice Thickness
- Do not Confuse Signal to Noise Ratio with Contrast to Noise Ratio!

### Some Theoretical Considerations

- Study Designs:
  - └ Blocked
  - └ Single Trial
- Predicting Responses
- Sources of Variance
- Resolution Limits:
  - └ Temporal
  - └ Spatial

### Blocked vs. Single Trial

**Typical Blocked Design**

25 sec  
Lights On  
Lights Off  
time

**Typical Single Trial Design**

A B C B A C Trial Type

Center for Cognitive Neuroscience

### Blocked Experiments

Stimulus  
Expected Response

- Responses are Fairly Large
- Data are Easy to Analyze
- With Long Blocks, Time course can be Ignored
- All trials within a block are treated as Identical

Center for Cognitive Neuroscience

### Single Trial Designs

A B C B A C Stimulus  
Expected Response

- Responses are small
- Useful contrast/noise is low
- Data are more Challenging to Analyze
- Exact Time course is Modeled or a Dependent Variable
- Suitable for Randomized Stimulus Designs

Center for Cognitive Neuroscience

### AFNI 3dToutcount

```
3dToutcount -automask KidAsImg.hdr | ldplot -stdin
```

Center for Cognitive Neuroscience

### P2P

Peak-to peak  
Std Dev  
Max  
Min

Center for Cognitive Neuroscience

### What really happened?

Center for Cognitive Neuroscience

- 
- HP filtering consequences
  - FT of scanner drift
  - Anti-alias at  $2^*tr$
  - linear regression for DCT gives fit values for the "nuisance" variables. Maps?
  - Gaussian smoother is a low pass filter
  - De-Meaning Figure